

Hive™ AUTOMATION PLATFORM

- *No downtime, automation runs continuously*
- *Load plates and consumables during runs*
- *Compact, versatile, and expandable*


The Hive automation platform automates laboratory processes and is fully customizable using standard products to match your workflow, storage and throughput requirements. The Hive footprint is similar to most laboratory benches so uses less of your precious lab space.

The unique design of the Hive allows use of the entire workspace, including above- and below-deck spaces. It facilitates compact system integration of BioNex instruments with third-party devices without compromising performance. Most devices face the front of the platform for easy offline use and serviceability. Not limited to a specific size, the Hive can be expanded to integrate multiple or larger devices.

Adding simultaneous-access storage allows plates and consumables to be loaded and unloaded during runs, eliminating traditional deck resetting time. One side of the storage module is dedicated to the Hive robotic arm while the other can be loaded by an operator or with BeeSmart Karts™, which are manually loaded at sample and consumable labs/stores remote from the Hive. Loaded Karts are driven to the Hive and docked for material transfer. The storage module is fully customizable for any combination of plate types and consumables at any capacity.



Hive Automation Platform

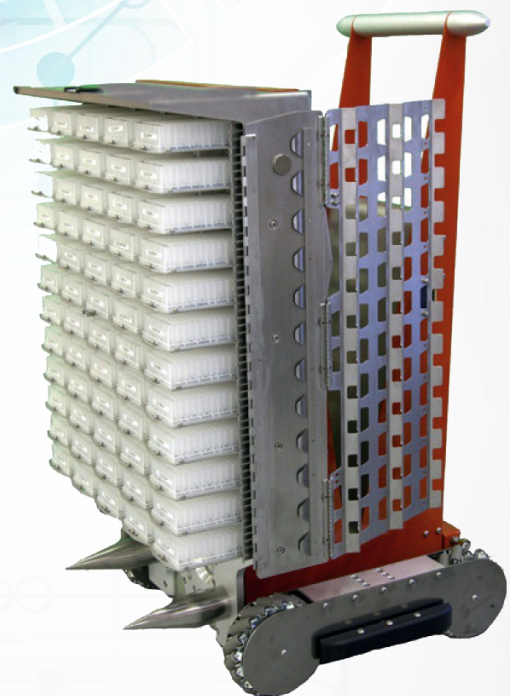

BeeSmart™ Kart



Hive Robot™


- 4-axis Hive Robot enables full utilization of Hive workspace
- Easy integration of third party instruments
- Robotic arm extends beyond Hive for additional workspace
- Fast, efficient pick and place across entire robot space
- Portrait and landscape labware handling
- Customizable gripper handles labware reliably
- Integrated imaging barcode reader with picture identification for reliable plate reading and tracking

Key Features and Benefits

- Increased workflow productivity
- Simultaneous-access storage
- Restock fresh and unload spent consumables during operation. Docking BeeSmart Kart
- Save lab space with compact Hive integrations
- Automate nearly any application by integrating a wide array of devices
- Forward facing instruments enable easy offline use and service
- Customizable random-access storage for consumables

General Specifications

Dimensions	System Design Specific
Power Source	100–240 VAC, 50/60 Hz, 2 A
Integrated UPS	Runs system for 45 minutes
Software	BioNex Synapsis™ Control Software

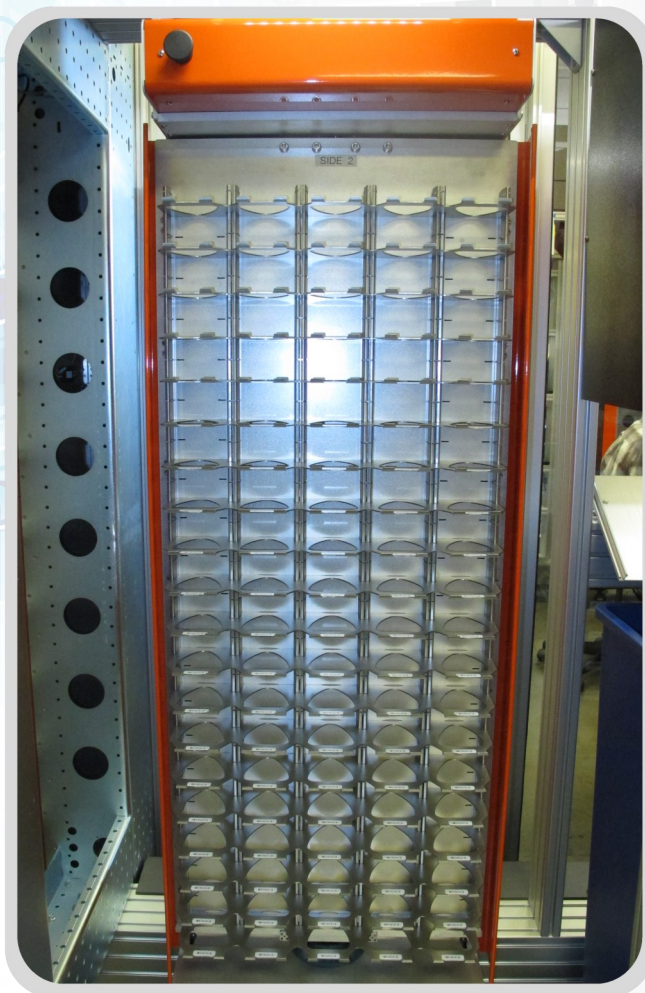

Simultaneous-access Plate Storage Module

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

LIT-50006 Rev B Effective 2017-01